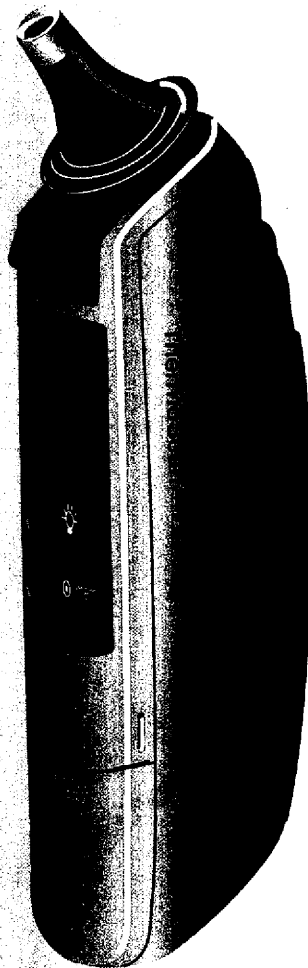
BRAUN
PRO 3000


Type 6014

ThermoScan®

English	4
Français	17
Español	31

Internet:
www.braun.com

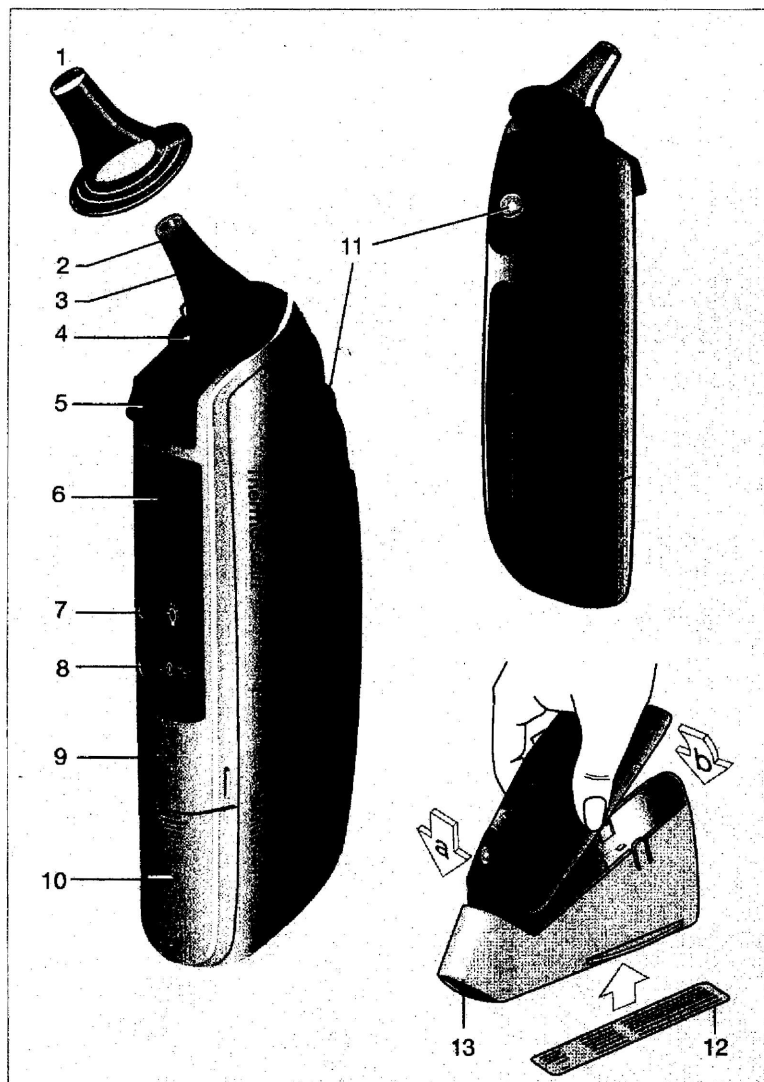6-014-101/01/III-00
USA/CDN/Mex.
Printed in Germany
Imprimé en Allemagne
Impreso en Alemania

English

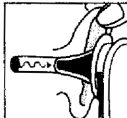
The Braun ThermoScan thermometer has been carefully developed for accurate, safe and fast human body temperature measurements in the ear. The shape of the thermometer prevents it from being inserted too far into the ear canal to damage the eardrum.

However, as with any thermometer, proper technique is critical to obtaining accurate temperatures. Therefore, please read all instructions carefully and thoroughly before using this product.

Important

- The probe window must be kept clean, dry, and undamaged at all times to ensure accurate readings. To protect the probe window, always keep the thermometer in the storage cover while transporting or when not in use.
- This thermometer must only be used with genuine Braun ThermoScan probe covers.
- Braun ThermoScan probe covers are for single use only.
- To avoid inaccurate readings, always attach a new, clean probe cover for each temperature reading.
- Always store the thermometer with a probe cover in place.

How does Braun ThermoScan work?

The Braun ThermoScan thermometer measures the infrared heat generated by the eardrum and surrounding tissue. To help ensure accuracy, the thermometer "scans" by taking 8 measurements in just one second and displaying the highest temperature.

The displayed ThermoScan temperature is the actual measured ear canal temperature plus a mathematical adjustment to approximate the familiar oral range. However, this is not necessarily the same as an oral temperature measured at the same time.

Why measure in the ear?

Clinical studies have shown that the ear is an excellent site for temperature measurement because temperatures taken in the ear reflect the body's core temperature. Body temperature is regulated by the hypothalamus, which shares the same blood supply as the tympanic membrane. Changes in core body temperature, are usually seen sooner at the tympanic membrane than at other sites, such as the rectum, mouth or under the arm.


Advantages of taking temperatures at the ear versus traditional sites:

- Axillary temperature readings only reflect skin temperature which may not indicate the internal body temperature.
- Rectal temperatures often lag significantly behind internal body temperature changes, especially at times of rapidly changing temperatures.
- Oral temperatures are often influenced by eating, drinking, thermometer placement, breathing through the mouth, or the inability of the person to close their mouth completely.

Body temperature

Normal body temperature is a range. The following table¹ shows that ranges of normal also vary by site. Therefore, readings from different sites, even if taken at the same time, should not be directly compared.

Axillary:	94.5°F - 99.1°F	34.7°C - 37.3°C
Oral:	95.9°F - 99.5°F	35.5°C - 37.5°C
Rectal:	97.9°F - 100.4°F	36.6°C - 38.0°C
ThermoScan ² :	96.4°F - 100.4°F	35.8°C - 38.0°C

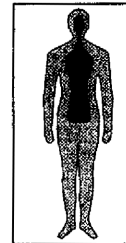Also, a person's normal temperature range tends to decrease with age. The following table shows normal ThermoScan ranges by age.

Normal ThermoScan temperature ranges²

0 - 2 years	97.5°F - 100.4°F	36.4°C - 38.0°C
3 - 10 years	97.0°F - 100.0°F	36.1°C - 37.8°C
11 - 65 years	96.6°F - 99.7°F	35.9°C - 37.6°C
> 65 years	96.4°F - 99.5°F	35.8°C - 37.5°C

The range of normal varies from person to person and can be influenced by many factors such as time of day, level of activity, medications and gender.

1. Chamberlain, J.M., Terndrup, T.E., New Light on Ear Thermometer Readings, *Contemporary Pediatrics*, March 1994.

2. Chamberlain, J.M., et al., Determination of Normal Ear Temperature with an Infrared Emission Detection Thermometer, *Annals of Emergency Medicine*, January 1995.

Product description

1. Probe cover
2. Probe window
3. Probe
4. Probe cover detector
5. Probe cover ejector
6. Display (LCD)
7. LCD light button ✱
8. Ⓞ/mem button (On/memory function)
9. Battery door lock
10. Battery door
11. Activation button
12. Label
13. Storage cover

Package components

- Braun ThermoScan thermometer
- Use and Care manual
- Quick Reference Guide
- Probe covers (20, plus one on the thermometer)
- Owner registration/warranty card
- Additional items may be included as noted on outer packaging
- Label
The first time you use the thermometer, please make sure to apply the special label included in the package, in the language of your choice (Item 12 on page 3).
- Infrared temperature reference available separately

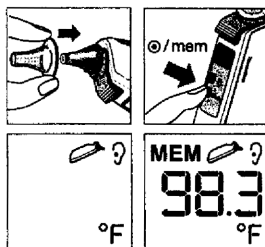
How to use Braun ThermoScan

1. **Always use a new, clean probe cover for each measurement to help ensure an accurate reading.**

Attach a new, clean probe cover or, if already attached, press the Ⓞ/mem button.

The LCD (liquid crystal display) is activated, showing all segments.

When the ready symbol  is displayed the thermometer is ready for use.



To help ensure an accurate temperature reading, use the following technique:
If you are right handed, hold the thermometer in the right hand and take the temperature in the right ear. If you are left handed, hold it in the left hand and use the left ear.

2. Perform an ear tug to straighten the ear canal. This gives the thermometer a clear view of the eardrum.

Children under 1 year:
Pull the ear **straight back**.

Children aged 1 year to adult:
Pull the ear **up and back**.

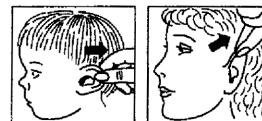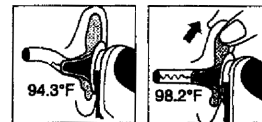
An ear tug is best performed by using your free hand to grasp the outer edge of the top half of the ear.


3. While tugging the ear, fit the probe snugly into the ear canal as far as possible and press the activation button. Release it when you hear a beep. This is the Temp Beep that confirms the end of measurement.

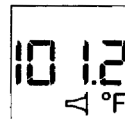
4. Remove the thermometer from the ear canal. The LCD displays the temperature measured and the probe cover symbol.

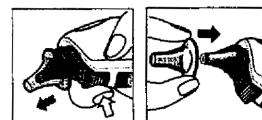5. **Replace the probe cover after each measurement:**
Press the ejector button and put on a new, clean probe cover.

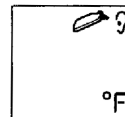Only then, the probe cover symbol disappears from the display. The thermometer is now ready for use again.

6. If no button is pressed or no new probe cover is attached within 30 seconds, the thermometer will turn off automatically.

Note: In the following situations, it is recommended that you take three temperatures in the same ear. If they differ, use the highest reading.

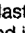
- Infants in the first 90 days of life.
- Children under three years of age who have a condition such as a compromised immune system and for whom the presence or absence of fever is critical.
- When you are first learning to use the ear thermometer until you are comfortable with the technique and are obtaining consistent readings.

Important: As with any type of thermometer, slight temperature variations (± 0.3 - 0.5°F / ± 0.2 - 0.3°C) can occur, if consecutive measurements are taken.

Temperature taking hints

- The right ear reading may differ from the reading taken at the left ear. Therefore, always take the temperature in the same ear.
- The ear must be free from obstructions or excess earwax build-up to take an accurate reading.
- External factors may influence ear temperatures, including when an individual has:
 - been lying on one ear or the other
 - had their ears covered
 - been exposed to very hot or very cold temperatures, or
 - been recently swimming or bathing.In these cases, remove the individual from the situation and wait 20 minutes prior to taking a temperature.


Memory mode

The last temperature taken before the thermometer powers down is stored in memory. To enter the memory mode, press the /mem button.

The last stored temperature is displayed along with the MEM symbol.


Even in memory mode, a new temperature can be taken provided that the ready symbol  is shown.

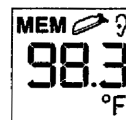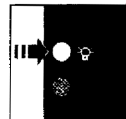To quit the memory mode, press the /mem button again.

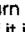
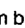
LCD light

for easy nighttime reading

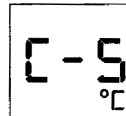
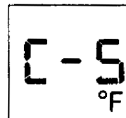
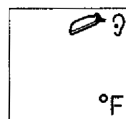
Changing the temperature scale

This thermometer is shipped with the Fahrenheit temperature scale activated. If you wish to switch to Celsius (°C), proceed as follows:

- Turn on the thermometer.
(If it is already turned on, make sure it is not in memory mode.)
- Press the /mem button and keep it pressed. Then press and release the activation button to switch over to the «change temperature scale» mode.
- By pressing the activation button again, the Celsius scale is activated, «°C» is displayed on the LCD. Each further pressing of the activation button toggles the scale between °C and °F.
- Release the /mem button to return to the ready symbol .

The temperature stored in memory will automatically be converted to the selected scale when displayed on the LCD.

Care and cleaning

Braun ThermoScan probe covers are intended for single use only! To ensure accurate readings, always use Braun ThermoScan disposable probe covers with the Pro 3000 thermometer. Damaged, perforated, or soiled probe covers can lead to inaccurate temperature readings. Additionally, the accuracy of your temperature readings can be affected by damage to the probe window, or the presence of dirt and earwax on the probe window.

Probe Window

To assure a high level of accuracy, it is very important to keep this window clean, dry and undamaged. Fingerprints, earwax, dust and other soiling compounds reduce transparency of the window and result in lower temperature readings. If the thermometer is accidentally placed in the ear without a probe cover installed, the window must be cleaned immediately. Additionally, the window and/or probe tip area should never be modified, changed, or adjusted. These changes will affect calibration and accuracy of the thermometer.



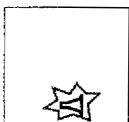
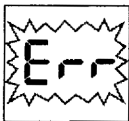
To clean the window, gently wipe its surface with a cotton swab slightly moistened with alcohol and immediately wipe dry with a clean cotton swab. After cleaning, allow at least 5 minutes drying time before taking temperatures. Make sure window is clean and dry. Avoid touching the window except when cleaning is required. If the window is damaged, the thermometer must be returned to Braun ThermoScan for service.




Note: Do not use any chemical other than alcohol to clean the probe window.

Thermometer


Use a soft cloth slightly moistened with alcohol to clean the thermometer display and exterior. Do not use abrasive cleaners. Never submerge the thermometer in water or any other liquid. Do not autoclave or immerse. Do not use ammonia or ammonium chloride based agents as they could damage the exterior of the instrument. If you currently use a specific cleaning agent or disinfectant, we recommend that you examine its chemical ingredients prior to use on the thermometer.

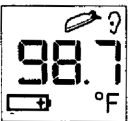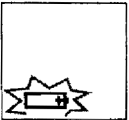Store thermometer and probe covers in a dry location free from dust and contamination and away from direct sunlight. The ambient temperature at the storage location should remain fairly constant and within the range of 50°F to 104°F (10°C - 40°C).

Troubleshooting		
Error message	Situation	Solution
 (J)	No probe cover is attached and activation button was pressed while probe was in the ear.	Make sure probe window is clean; refer to section «Care and cleaning» if necessary. Attach new, clean probe cover to stop error beeps.
	Ambient temperature is not within the allowed operating range of 50°F to 104°F (10°C-40°C) or changing too rapidly.	Allow the thermometer to remain in a room for 30 minutes where the temperature is between 50°F to 104°F (10°C-40°C).
	Temperature taken is not within displayed temperature range (68°F-108°F / 20°C-42.2°C).	Make sure new, clean probe cover is attached and thermometer is properly inserted. To clear the error message, attach a new, clean probe cover or press the ©/mem button.
	HI = too high LO = too low	Then, take a new temperature.
	Battery is low, but thermometer will still operate correctly. (No light function)	Insert new batteries
	Battery is too low to take temperature readings.	Insert new batteries
	System error	Wait 30 seconds until thermometer powers down, then turn on again. ... reset the thermometer by removing the batteries and putting them back in. ... call Braun ThermoScan Customer Service at 1-800-327-7226.
	• If error persists	
	• If error still persists	

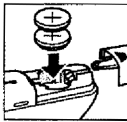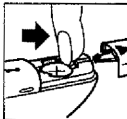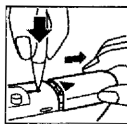Replacing the batteries

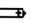The thermometer is supplied with two 3 V lithium batteries (CR 2032/DL 2032). Insert new batteries when the battery symbol  appears on the display. Using the tip of a ball-point pen, press the battery door lock to open the battery compartment.

Using your finger, press the battery release.

Remove the batteries and replace with new batteries, making sure the poles are in the right direction. Slide battery door back in until it snaps in place.

Do not dispose of empty batteries in regular waste. Take them to special collection sites.

Special situations

Clinical studies have shown the Braun ThermoScan thermometer can obtain accurate temperatures on persons of all ages. However, there are certain situations when an ear thermometer should not be used. These include, but may not be limited to, the following situations:

- Do not use an ear thermometer if there is blood or drainage in the external ear canal.
- An ear thermometer should not be used on any patient's ear which exhibits symptoms of an acute or chronic inflammatory condition of the external ear canal.
- Common situations like the presence of moderate amounts of ear wax in the ear canal, otitis media and tympanostomy tubes do not significantly impact temperature readings. However, complete ear canal occlusion due to ear wax can result in lower temperature readings.
- For persons wearing hearing aids or ear plugs, remove the device and wait 20 minutes prior to taking a temperature.
- Use the untreated ear if prescription ear drops or other ear medications have been placed in the ear canal.
- Patients who have deformities of the face and/or ear may not be able to have a temperature taken with an ear thermometer.

Fever facts

Many persons may not have an elevated temperature even when they are ill. These include, but are not limited to infants under 90 days of age, persons taking steroids, antibiotics or antipyretics (acetaminophen, ibuprofen, aspirin), persons with compromised immune systems, including the elderly or persons with some chronic illnesses.

Fever is described as an elevation of body temperature over an individual's «normal» temperature.

An elevated temperature or fever is often viewed as a danger sign. In fact, a fever can be very beneficial, and helps our immune system work more effectively. A fever should be evaluated in the light of other physical symptoms. With the exception of newborn infants, the presence of fever, without any other symptoms of illness, or in a child who is behaving normally may not be cause for concern.

Some people, like the elderly, infants under 90 days of age, those with compromised immune systems, or persons who take steroids, for example are often unable to build a response to illness or environmental conditions. These individuals may not have a fever when they are ill, or their elevated temperature may be lower than expected for the severity of their illness. Other medications such as anti-inflammatories and some analgesics may also mask fever.

The presence or absence of fever should rarely be used as the only measurement of the significance of illness.

Antipyretics, like acetaminophen or ibuprofen, are usually recommended to relieve the aches and associated symptoms of fever, not the fever itself. Febrile seizures, or convulsions, which usually occur in children six months to six years of age, are thought to occur not because a fever is present, but because of the rapid rise of the child's temperature.

Product specifications

Displayed temperature range	68 °F - 108 °F (20 °C - 42.2 °C)
Operating ambient temperature range	50 °F - 104 °F (10 °C - 40 °C)
Display resolution	0.1 °F or °C

Accuracy for patient temperature range

	<u>Maximum Laboratory Error</u>	
96.8 °F - 102.2 °F (36 °C - 39 °C):	±0.4 °F	±0.2 °C
outside this range:	±0.5 °F	±0.3 °C

Long term storage ranges

Temperature:	-4 °F to 122 °F (-20 °C to 50 °C)
Humidity:	95% non-condensing

Shock: Withstands one meter drop

Measurement time	1 second
Minimum time between readings	2 seconds
Automatic «power down»	30 seconds

Battery life	5,000 temperatures
Battery type	2 x CR 2032/DL 2032
Thermometer dimensions	6.4" x 1.7" x 1.2" (16.3 cm x 4.4 cm x 3.1 cm)
Thermometer weight (without batteries)	3.5 ounces (100 g)

This infrared thermometer meets requirements established in ASTM Standard E 1965-98. Full responsibility for the conformance of the product to the standard is assumed by Braun GmbH, 61476 Kronberg, Germany.

ASTM laboratory accuracy requirements in the display range of 96.8 °F to 102.2 °F (36 °C to 39 °C) for infrared thermometers is ±0.4 °F (±0.2 °C), whereas for mercury-in-glass and electronic thermometers, the requirement per ASTM Standards E667-86 and E1112-86 is ±0.2 °F (±0.1 °C).

A clinical summary is available upon request.

The Braun ThermoScan thermometer has been proven to be safe. To ensure the best results, please take the time to read this manual completely and keep it handy for future reference. In the unlikely event you experience any difficulty using the Braun ThermoScan thermometer, simply call us toll free at: 1 (800) 327-7226.

This equipment has been tested for potential electromagnetic interference in accordance to IEC 601-1-2. These tests confirm that the emissions generated by the Braun ThermoScan Pro 3000 were well below the required limits. Thus, interference with other equipment is not to be expected.

TYPE BF
EQUIPMENT

Attention, consult
ACCOMPANYING
DOCUMENTS

MEDICAL ELECTRICAL EQUIPMENT
CLASSIFIED BY UNDERWRITERS LABORATORIES INC.®
WITH RESPECT TO ELECTRIC SHOCK, FIRE AND
MECHANICAL HAZARDS ONLY
IN ACCORDANCE WITH UL 2601-1 / CAN CSA C22.2
No. 601.1 10HK

Internally Powered Equipment
Continuous Operation
Not Protected against Ingress of Water
Patents Pending

Service

The Braun ThermoScan thermometer has a limited three year warranty. Should service be required during or after the warranty period, you must call Braun ThermoScan's Customer Service Department at 1-800-327-7226 to obtain a return authorization number and shipping address. Repackage the thermometer carefully in its original packaging or securely packed to avoid damage during shipping. Include the original sales slip indicating the date of purchase, a note describing the problem, your authorization number, and your return address. Send the Braun ThermoScan thermometer prepaid and insured.

Note: Shipping address varies by country.

Answers to any questions regarding operation of the Braun ThermoScan thermometer may be obtained by calling us toll free at:

(800) EAR SCAN
(800) 327-7226
Monday-Friday
6:30 A.M. to 5:00 P.M. Pacific Time
www.braun.com

Limited three year warranty

For Model Pro 3000

The Pro 3000 instrument is warranted by Braun to be free from defects in material and workmanship under normal use and service for a period of three years from the date of delivery to the first professional user who purchases the instrument.

This warranty does not cover batteries, damage to the probe window, or damage to the instrument caused by misuse, negligence or accident, and extends only to the first purchaser of the product. Additionally, this warranty becomes void if the thermometer is operated with anything other than genuine Braun ThermoScan® Probe Covers.

Thermometer service or replacement will be carried out at Braun's expense subject to the terms of this warranty. Loss or damage in return shipment to Braun shall be at purchaser's risk. Braun reserves the right to replace a defective product with the most comparable product currently available. This product has a five year minimum expected lifetime.

IN NO EVENT SHALL BRAUN BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH THE PURCHASE OR USE OF ANY PRODUCT. This warranty shall not apply to, and Braun shall not be responsible for, any loss arising in connection with the purchase or use of any product which has been serviced by anyone other than a Braun authorized service representative or altered in any way so as, in Braun's judgment, to affect its stability or reliability, or which has been subject to misuse or negligence or accident, or which has had the serial number or lot number altered, effaced or removed, or which has been used other than in accordance with the instructions furnished by Braun. IN NO EVENT SHALL BRAUN'S LIABILITY EXCEED THE ORIGINAL COST OF THE PRODUCT.

This warranty is in lieu of all other warranties, express or implied, and of all other obligations or liabilities on Braun's part, and Braun neither assumes nor authorizes any representative or other person to assume for it any other liability in connection with the sale of its products.

TO THE EXTENT PERMITTED BY LAW, BRAUN DISCLAIMS ALL STATUTORY AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO, WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

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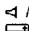
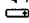
(USA)


Error Messages

Hi	temperature too high
Lo	temperature too low
Err	ambient temperature too low / too high / changing too fast
all symbols flashing	system error
 	old / no probe cover attached
	battery low

ThermoScan ear thermometer For use in ear canal only
Normal ThermoScan body temperature range
 35.8 - 38.0 °C 96.4 - 100.4 °F
Key Operating Points
 For detailed instructions refer to Use and Care Manual 
 • Keep probe window clean at all times
 • Always attach a new clean probe cover
 • Perform an ear tug
 • Fit probe snugly into ear canal
 • Press and hold activation button until you hear Temp Beep

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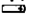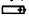
Brünn

Body temperature

- Normal body temperature is a range. The following table shows that this normal range also varies by age. Therefore, readings from different sites should not be directly compared.

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Normal ThermoScan temperature ranges?

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1. Always use a new, clean probe cover for each measurement to help ensure an accurate reading.
Attach a new, clean probe cover or, if already attached, press the  /mem button.

2. Pull an ear tug to straighten the ear canal. This gives the thermometer a clear view of the eardrum.

Children under 1 year:
Pull the ear straight back.
Children aged 1 year to adult:
Pull the ear up and back.

3. While tugging the ear, fit the probe snugly into the ear canal as far as possible and press the activation button. Release it when you hear the Temp Beep.

6-014-027000-09